


Pumpsure



INFUSION PUMP

Easy To Learn Easy To Use

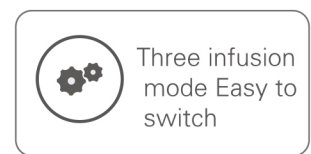

The pumpsure 5000 series syringe pump has a friendly interface. The operation is easy and fast. There are dual system architecture and double pressure monitoring system which will guarantee the safety of the whole infusion process. It supports multiple infusion modes including rate mode, time mode and weight mode.

The pump can be used in ICU, anesthesiology, pediatrics, oncology, general unit and emergency department.



SPECIFICATIONS

Infusion set brand	Support specified 4 brands, 4 custom brands	<div style="border: 1px solid black; padding: 5px; display: inline-block;">IPX4</div> Splashproof
Infusion set size	20 d/ml, 60 d/ml	
Infusion rate	0.1~1200mL/h	
Volume range	0.1-9999ml	
Air detection	High, medium, low	
K.V.O	off, 0.1ml/h-5ml/h	
Bolus Rate	0.1-999ml/h	
Nearly empty time	0-15min	
Accuracy	±5% (Measured at a 25 mL/h infusion rate)	
Infusion mode	Time mode, rate mode, weight mode	
Occlusion level and range	High, medium, low, range: 50~1300mmHg	
Weight	1.8kg (with pole clamp)	
Dimension	290*120*125(mm)	
Data storage time	Setting & infusion history can be kept for at least 10 years after shutdown	
Battery type	DC 11.1V recharge battery pack, after fully charged, continuous infusion is more than 7 hours at a 25 mL/h infusion rate	
Power input	AC 100-240 VAC	
Alarms	Battery empty; Motor rate error; Downstream occlusion; Upstream occlusion; Air Alarm; Drop speed abnormal; Near empty; Empty; Low battery; Reminder alarm	

The above specification is only suitable for 5000BP. As for other models, please refer to the user manual.